







MODEL & SPECIFICATION

Disposable non-absorbable ligation clips

MODEL	SPECIFICATION	CLIPS/CARTRIDGE	CARTRIDGES/BOX	LENGTH	CLOSING RANGE
 Type I	Small	2 / 4 / 6	14	9.5mm	1mm-3mm
 Type II	Medium	2 / 4 / 6	14	13mm	2mm-7mm
 Type III	Large	2 / 4 / 6	14	17.5mm	3mm-10mm

Clip applicator/Clip remover

MODEL OF CLIP	MODEL OF CLIP APPLIER	MODEL OF CLIP REMOVER
 Type I	KSQ1, KSK1	KCQ1, KCK1
 Type II	KSQ2, KSK2	KCQ1, KCQ2, KCK1, KCK2
 Type III	KSQ3, KSK3	KCQ2, KCK2

Disposable continuous clip applicator

MODEL	SPECIFICATION	MODEL OF CLIP	NUMBER OF CLIPS PER UNIT
2711005	Small	Type I	5
2711006	Small	Type I	6
2711008	Small	Type I	8
2711010	Small	Type I	10
2711012	Small	Type I	12
2711013	Small	Type I	13
2711015	Small	Type I	15



DISPOSABLE NON-ABSORBABLE LIGATION CLIPS

DISPOSABLE CONTINUOUS CLIP APPLIER

Condiner Medical Technology Co., Ltd.

No.117 Xueye Road, Xinbei Zone, Changzhou, Jiangsu, China, 213125
Tel : 86-519-68585073
Fax : 86-519-68585080
Web: www.condiner.com.cn



PRODUCT PERFORMANCE


- Disposable non-absorbable ligation clips are applied to close the tubular tissue structures like cystic duct, arterial vascular or venous. Its principle is to achieve the effect of ligation or closure of the vascular and cystic duct tissue by using this product.

SCOPE OF APPLICATION

- Applied to close the tubular tissue structures like cystic duct, arterial vascular or venous. Compatible with Condiner clip applicator model KSQ1,KSQ2,KSQ3,KSK1,KSK2,KSK3.

PRODUCT ADVANTAGES

- Meet all clinical demands with large, medium and small models
- Suitable for all kinds of surgeries, easy, fast, safe and effective
- High polymer material with good biocompatibility
- Non-slip tooth design, effectively prevent sliding
- Design of safety lock
- Accurate clipping confirmation with tactile feedback
- Pressure adaptive design at the end and flexible hinge ensures closure without cutting tissues
- Radiolucent
- Revokable ligation

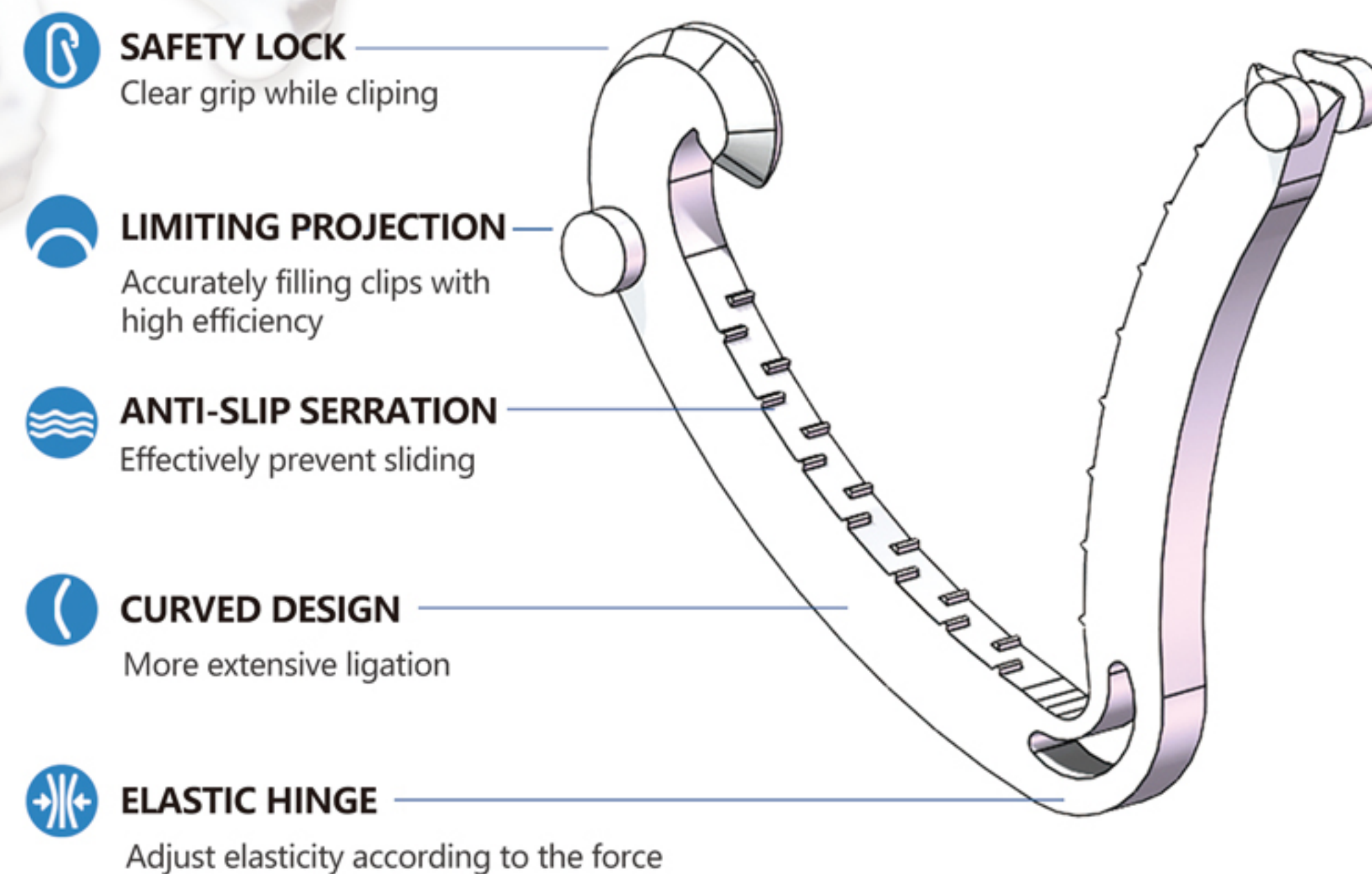


PRODUCT PERFORMANCE

Disposable continuous clip applicator is sterilized medical device, which consists of ligation clips and continuous clip applicator. While pressing the compression bar, it will push out the clips automatically. When the rotating head is turned, the applicator head will also rotate accordingly.

SCOPE OF APPLICATION

Disposable continuous clip applicator is used to ligate vessels or tubular tissues during endoscopic surgery.

PRODUCT ADVANTAGES

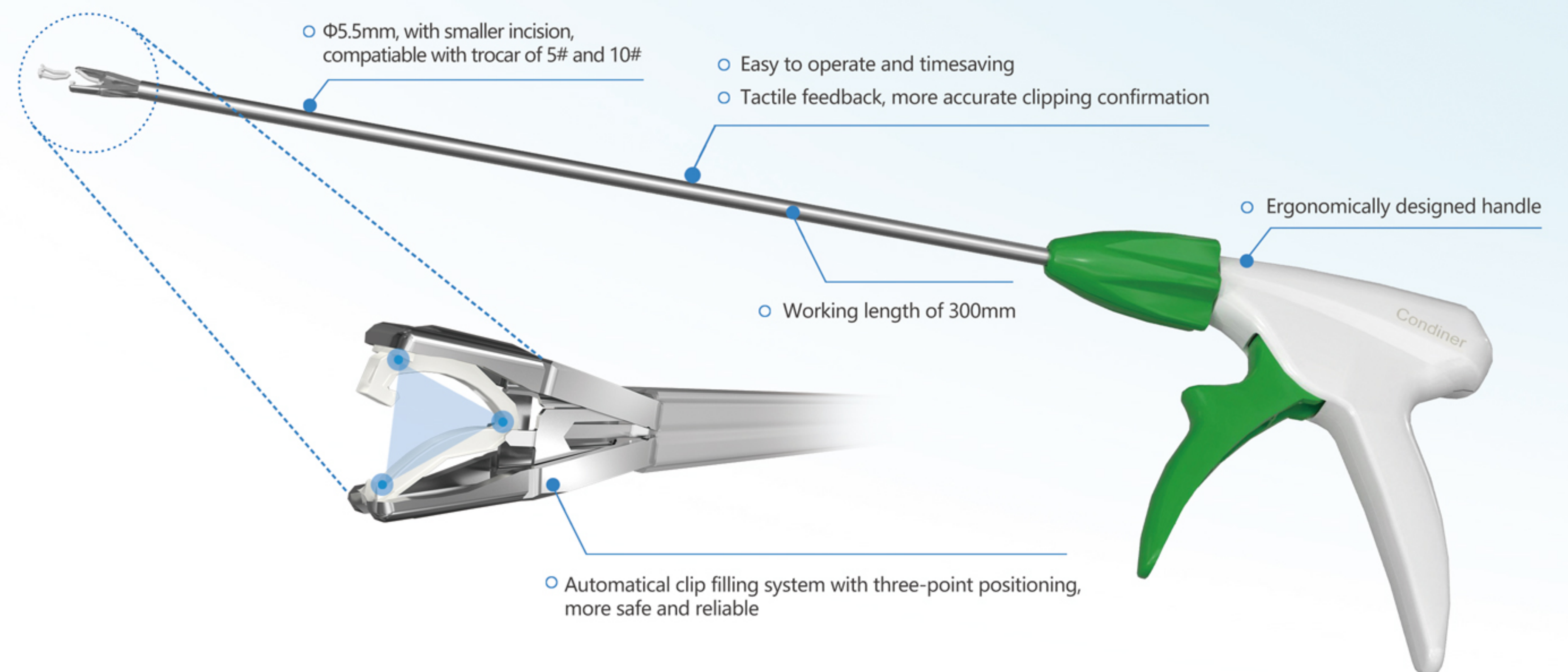

LAST PIECE REMIND

When last clip is ligated, the orange marker will occur to remind user to stop clipping